
BERDSEN

Smart Oral Irrigator ClearJet X3

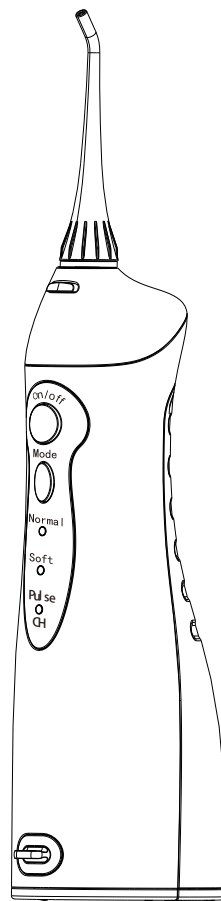EN

OPERATING MANUALS

KEEP THESE INSTRUCTION FOR FUTURE REFERENCE.

IMPORTANT SAFEGUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

READ ALL INSTRUCTIONS BEFORE USING

DANGERS!

- Toddlers, infants, people who are unable to operate should not allowed use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into tub or sink.
- Do not reach for the product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.


WARNINGS!

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery;
- The battery is to be disposed of safely.
- Regarding how to remove the battery, please to see "Disposal" in **page 7**.
- Do not fill the tank with water warmer than 40 °C.

- If the external flexible cable or cord of this transformer is damaged, it shall be replaced by the manufacturer or his service agent or a similar qualified person in order to avoid a Hazard.
- Keep the cord away from heated surfaces.
- Please do not use it if you lost the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be become lodged in the throat.
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums. Doing so may injure your gums.

SAVE THESE INSTRUCTIONS

CONFIGURATION

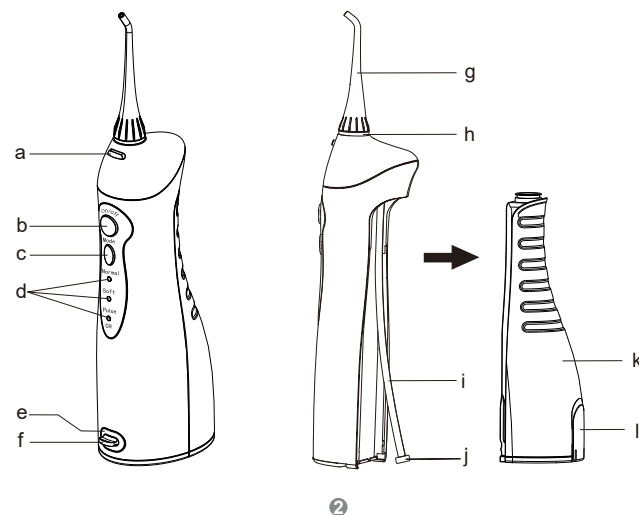- a. Nozzle release button
- b. On/Off button
- c. Mode button
- d. Mode indicator
- e. Charging base
- f. Charging base cover

- g. Nozzle
- h. Identification ring
- i. Nozzle suction tube
- j. Filter
- k. Water tank
- l. Water tank cap

Note:

There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

HOW TO USE

Charging the device(Fig. 1)

1. Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).
2. Insert the charging adaptor securely into the power outlet.
3. The "CH" indicator blinks when the device is charging, and will turn off after the device is fully charged.
4. Recover the cover onto charger base again after charging finished.

Note:

If the "CH" indicator is blinking during or after an irrigating, that means the battery need to be recharged. An full recharging will take about 8 hours.

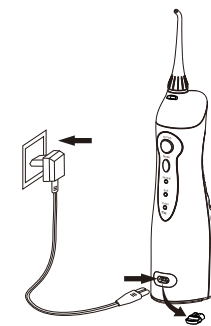
Fig.1

Use the Oral Irrigator

1. Attach nozzle onto the unit. (Fig. 2)
2. Open the water tank cap to fill the tank with water which temperature not exceed 40°C, then close the cap. (Fig. 3)
3. Or, Slide the water tank down off the unit, fill the tank with water which temperature not exceed 40°C, then attach the tank onto the unit. (Fig. 4)
4. Press "Mode" key to select the mode you desired. During irrigating, you can change the working mode. (Fig. 5)
Normal: for normal cleaning.
Soft: for gentle cleaning of sensitive gums.
Pulse: for cleaning and massaging the gums.
5. Hold the unit in an upright position with nozzle tip point toward your teeth.
6. Keep your mouth opened slightly so that the water can flow out.
7. Press "On/Off" key to start irrigating. (Fig. 6)

- 1). Directing the water jet flow at a 90-degree angle to the teeth and gums
- 2). Directing the water jet slowly along the teeth
- 3). Directing the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue
- 4). Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

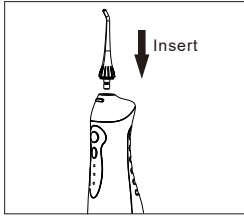
Fig. 2

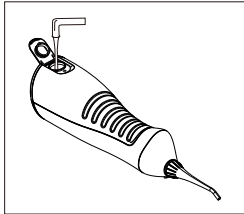
Fig. 3

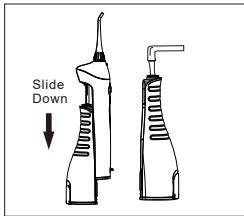
Fig. 4

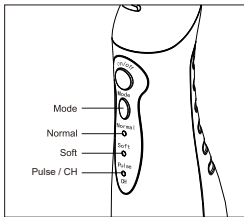
Fig. 5

8. Press "On/Off" key and then press the nozzle release button to replace other nozzle when you need. (Fig. 7)
9. Once irrigating complete, press "On/Off" key again to turn off the device.

After using the appliance (Fig. 8)

1. Open the water tank cap and dispose of the water in the tank.
2. Press "On/Off" key to turn on the unit, and drain the water inside the appliance.
3. Press "On/Off" key again to turn off the unit, and close the water tank cap.
4. Wipe the device with dry cloth.

Note:

To prevent bacteria propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for long time, be sure wipe and dry the appliance ASAP before store it.

MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 50°C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

1. Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

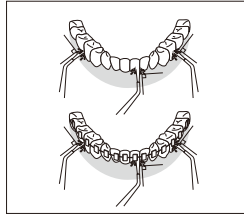
Fig. 6

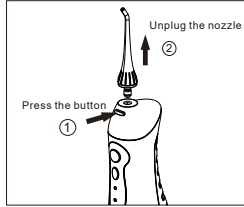
Fig. 7

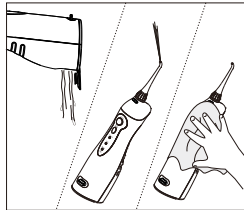
Fig. 8

Water tank cap

Wipe off any grime and dirt with a cloth.

Nozzle

1. Wash with water and wipe with soft cloth.
2. Replace the nozzle every 6 months.

Water suction tube

1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.

Tank

1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

TROUBLESHOOTING

Problem	Possible Cause	Solution
It does not work	You've just purchased the product or it has been left unused for more than 3 months.	Charge the unit for at least 8 hours.
	The battery has reached the end of its lift cycle.	Contact service center for support.
It only works for a few minutes even after charging	Charging time was insufficient.	Charge the unit for at least 8 hours.
	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
Water pressure is very low	There is dirt stuck in the filter	Clean the filter
	The tank is empty	Fill up tank
Water does not come out	The appliance is tilted too much while in use	Use with the appliance in a vertical state

DISPOSAL

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove 2pcs rubber feet on the bottom the handle. (Like ① on Fig. 9)
 2. Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (Like ② on Fig. 9)
 3. Slide the water tank down off the device and unscrew the two screws beside of the suction tube. (Fig. 10)
 4. Open the top cover and put the inner out of the unit, then unscrew all the screws of the inner. (Fig. 11)
 5. Cut off the battery line and pull out the batteries from the inner. (Fig. 12)
- If you have recycling questions, please contact your local waste management office.

BERDSEN

Made in China for:

Global Income sp. z o.o.
ul. Chwaszczyńska 135B
81-571 Gdynia, Poland

www.berdsen.pl.

IMPORTANT ENVIRONMENTAL INFORMATION ABOUT THIS PRODUCT

This symbol on the device or the package indicates that disposal of the device after its lifecycle could harm the environment. Do not dispose of the unit (or battery) as unsorted municipal waste; it should be taken to a specialized company for recycling. This device should be returned to a local recycling service. Respect the local environmental rules. If in doubt, please contact your local waste disposal authorities.

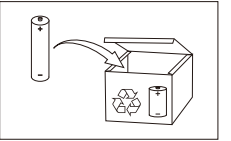
Fig. 9

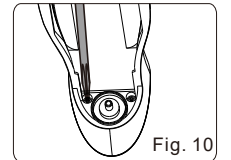
Fig. 10

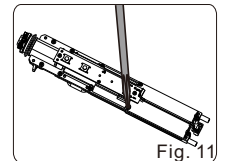
Fig. 11

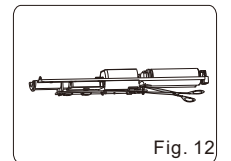
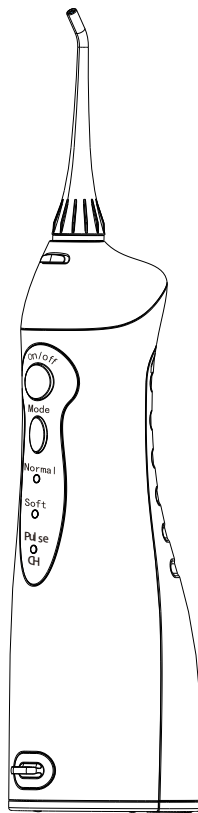
Fig. 12

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