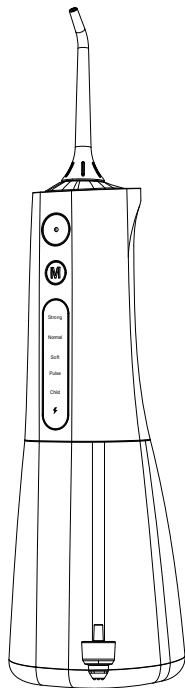
BERDSEN

BD-200

PORTABLE ORAL IRRIGATOR

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EN

Instruction manual

Read and follow the instruction before installation and use.
Keep the instruction for future reference.

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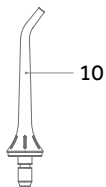
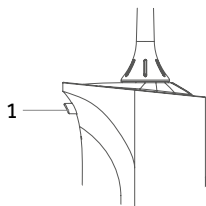
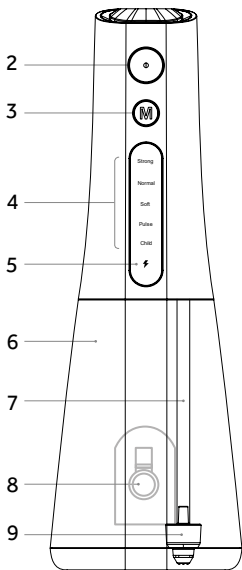
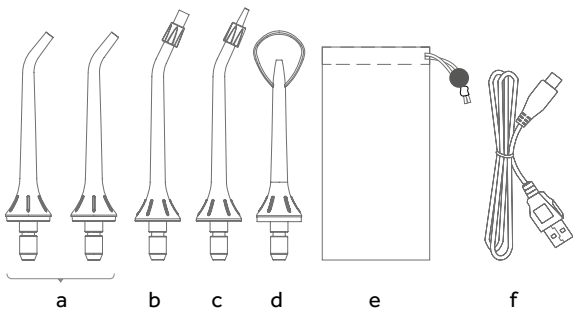
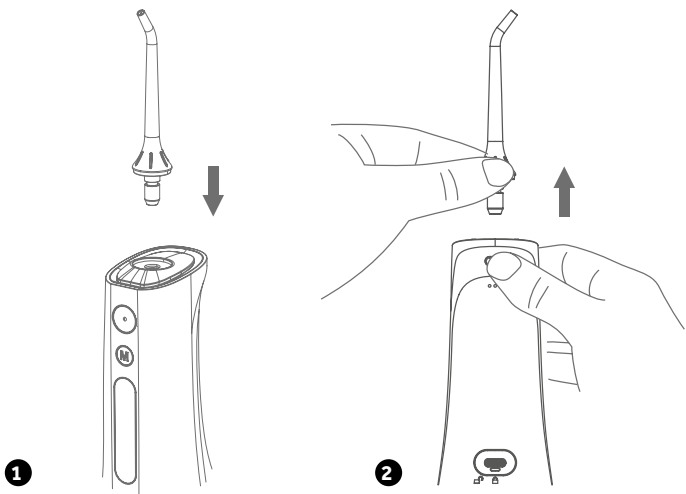
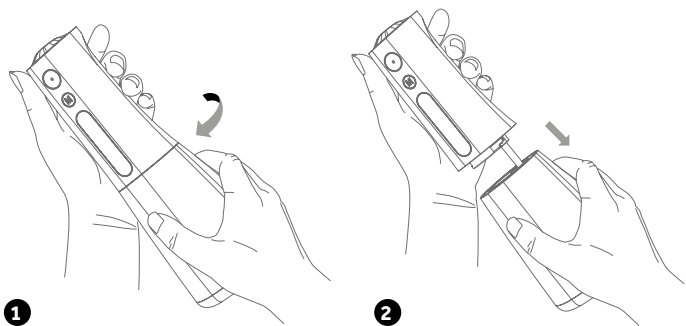
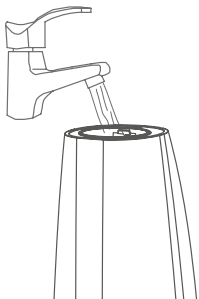
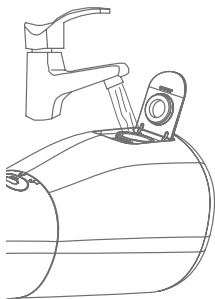
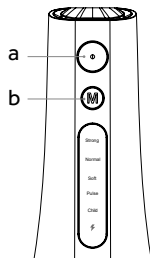
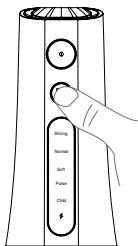
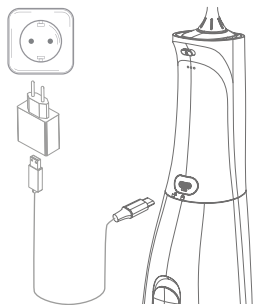
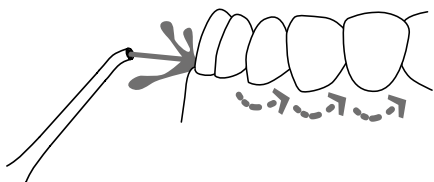
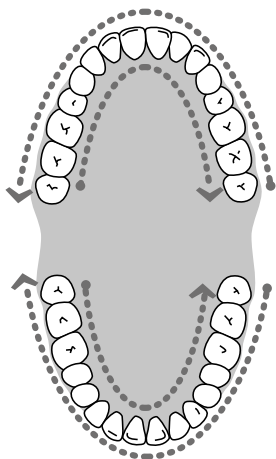
III**IV**

V**1****2****VI****VII****VIII****IX**

Thank you for choosing a Berdsen product. Please read this manual carefully before use and keep it for future reference.

SAFETY WARNINGS

- The appliance must not be used by children under 8 years of age. Prevent them from having access to the appliance.
- Before and after use, check the USB cable for damage and keep it away from heated surfaces.
- Do not use the device if the power cable is damaged.
- You can use your mobile phone charger power bank, computer USB port, and other power supplies to charge. However, never charge your device with a charger with an output voltage of more than 5V.
- Disconnect the power cable when you have finished using the device.
- Do not connect the power cable with wet hands.
- Do not place the device in water or other liquids.
- Do not place the device where there is a risk of it falling into water.
- Do not use the device while bathing.
- Do not touch the appliance if it has fallen into water. Disconnect the power cord from the mains socket immediately.
- Do not disassemble or repair the interior of this product yourself.
- Do not connect any foreign objects to the device.
- Do not touch the product's power supply while charging and expose it to adverse conditions.
- Do not place or store the product where it can fall or be pulled into drainage holes.
- Only fill the product with water.

WARNINGS

- Keep jet tips away from small children. The tips are small enough to be swallowed by them and there is a risk of choking.
- The appliance is not a toy.
- The appliance may only be used by children aged 8 years and over and by persons with limited physical, mental or sensory abilities, or persons without knowledge or experience of using this type of appliance, if supervised.
- If you have had oral or dental surgery within the last 2 months, consult your doctor and follow his advice before using this product.
- When using this oral cleansing product for the first time, slight bleeding of the gums may occur, even if the gums are healthy. This is due to the fact that the gums are being stimulated by the jet pressure for the first time. Usually the gums stop bleeding after 1-2 weeks of use. If the bleeding of the gums continues for more than 2 weeks, this may indicate pathological conditions, such as periodontitis. If this is the case, stop using the product and consult your dentist.
- Do not use the device if you experience pain in your teeth or gums. Consult your dentist.
- Do not put the tip too close to your teeth or gums, as you may damage the gums.
- Device for private home use, not intended for use in a dental practice.
- The appliance should only be connected by an adult.
- Plug the appliance into a grounded socket.
- Make sure the cord, plug and housing are whole and undamaged before plugging into an outlet.
- Ensure that the power cord does not come into contact with water.
- Do not extend or repair the power cord.
- Do not get the product wet or hot.
- Do not connect to more than one power supply.
- Do not use extension cords.

- The appliance will only be completely disconnected from the power supply when the plug is disconnected from the socket.
- Disconnect the plug from the socket before cleaning the product.
- Do not use the product if the USB cable or plug of the product is damaged, malfunctions, has fallen, been damaged or has fallen into any liquid.
- Only use the product with the power cable supplied.
- Do not use accessories other than those recommended by the manufacturer.
- Do not use the appliance without the jet tips.
- Do not direct the water jet directly under the tongue, into the eyes, nose or other sensitive areas. The pressure generated by this product can cause serious damage.
- Do not use the appliance for any purpose other than cleaning your teeth.
- Do not refill the water tank with mouthwash as it may cause damage.
- Fill the water tank with clean water.
- Do not fill the water tank with water above 40°C.

CONSTRUCTION OF THE DEVICE (FIG. I)

- | | |
|--------------------------|--------------------|
| 1. Jet tip unlock button | 6. Water tank |
| 2. On/off button | 7. Water hose |
| 3. Mode button | 8. Water tank plug |
| 4. Mode indicator | 9. Gravity ball |
| 5. Charging indicator | 10. Jet tip |

Note: Before first use, there may be residual fluids in the water tank and hose. This is the residue of the distilled water used during the factory inspection. This poses no risk to the health and safety of the user.

ACCESSORIES (FIG. II)

- a. standard tips
- b. orthodontic tip
- c. periodontal tip
- d. tongue cleaner
- e. pouch
- f. USB-C cable

BEFORE USE

INSTALLATION AND REMOVAL OF THE JET TIPS (FIG. III)

1. **Assemble:** Insert the desired jet tip into the mounting hole until you hear a clicking sound.
2. **Disassembly:** Press and hold the jet tip unlock button and at the same time pull the tip upwards to remove it.

Note: Do not press the jet tip unlock button during use to avoid accidental injury due to the tip falling out of the mounting hole under water pressure.

INSTALLATION AND REMOVAL OF THE WATER TANK (FIG. IV)

1. Grab the oral irrigator body with your left hand and the water tank with your right hand.
2. Turn the water tank clockwise.
3. Remove it by pulling downwards.

To reinstall the water tank, follow the above steps in reverse order.

FILLING THE WATER TANK (FIG. V)

The water tank can be filled as follows:

1. Through the top opening in the water tank (Fig. V.1) - to do this, dismantle the water tank according to the guidelines in this manual.
2. Through the opening at the back of the water tank (fig. V.2) - Open the water tank plug, fill the tank with water and then close the plug. In this case, the water tank must be fitted to the irrigator.

Note: Never use the appliance without a water tank or with an empty water tank. Only fill the irrigator with clean warm water under 40°. It is forbidden to add salt, toothpaste, mouthwash or other chemicals.

EMERGENCY ADDING OF WATER (FIG. VI)

When using the oral irrigator for the first time or after not using the oral irrigator for a long time, it may take a while for the jet tip to release water when the on/off button is pressed. If the water jet does not appear, remove the tip, add a little water to the mounting hole, reinstall the tip and then switch on the device.

MODE SELECTION (FIG. VII)

Press the on/off button to switch the irrigator on.

Then press the mode button until you select the desired operating mode.

Note: The oral irrigator is equipped with an additional protection mechanism. After two minutes of use, the device will switch off automatically. To start it up again, press the on/off button.

OPERATING MODES

Strong	Strong water pressure, suitable for deep cleaning and severe impaction.
Normal	Standard water pressure, high cleaning performance, suitable for daily teeth cleaning.
Soft	Low water pressure, suitable for sensitive teeth and new users.
Pulse	Strong and weak alternating water pressure, comfortable use, good for gum protection and daily cleaning.
Child	Gentle water pressure, suitable for children.

BATTERY CHARGE STATUS (FIG. VIII)

You can check the battery charge status by pressing the mode button for a few seconds.

The number of mode indicators that light up will indicate the battery charge level:

Number of lit mode indicators	Battery charge level
5	80-100%
4	60-80%
3	40-60%
2	20-40%
1	less than 20%

When the oral irrigator is in use and the battery needs charging, the charging indicator  will light up red for 15 seconds and then switch off.

CHARGING (FIG. IX)

Charge the oral irrigator for 6 hours before first use.

1. Switch the oral irrigator off using the on/off button.
2. Connect one end of the USB charging cable to the charging port in the oral irrigator and the other end to the 5V charger.
3. Plug the charger into a power socket.
4. During charging, the charging indicator  will light up red and when fully charged, the indicator will light up green.

Notes:

You can connect the USB cable to all adapters with max output up to 5V, e.g.: mobile phone charger, power bank, computer USB port and other power supplies.

When charging, place the device in a stable position.

Make sure the charging port is clean and dry.

USING THE ORAL IRRIGATOR (FIG. X)

1. Place the jet tip end in the mouth.
2. Press the on/off button to switch the device on.
3. Select the desired operating mode and aim the water jet at the teeth, the gap between the teeth or the gums.
4. Make sure that the water jet is perpendicular to the parts to be cleaned.
5. You can start cleaning from the molar part (back teeth) and continue towards the front teeth. Move the jet tip along the gum line, paying particular attention to cleaning the gaps between teeth and the area around the braces.
6. Press the on/off button to switch the device off.
 - For people with sensitive teeth, it is recommended to select a weaker water pressure.
 - When in use, hold the device in a position that allows water to be drawn from the water tank through the water hose, which is supposed to be facilitated by the gravity ball attached to the water hose. It is advisable to hold the device upright so that the jet end is above the water tank. Do not turn the device 180°.

Note: Keep your mouth slightly open to allow the water to flow freely.

7. Remove the water tank, pour out the remaining water in it, wipe off the moisture and place the device in a dry place.

CLEANING AND MAINTENANCE

Use only clean water and natural cleaning agents to clean the appliance. Do not use strong detergents or chemicals that may damage the appliance. Do not wash the appliance with water above 50°C. Keep the device away from heat sources and avoid direct contact with the sun.

MAIN UNIT

1. Disconnect the water tank and remove any residual water from the inside of the unit.
2. Clean the unit with a soft cloth and liquid soap and rinse it thoroughly and dry it with a dry cloth.
3. Do not immerse the appliance in water.
4. Do not use alcohol, petrol, thinner etc. to wipe down the product.

WATER TANK

1. Wash the water tank and dry it with a dry cloth.
2. When the device is not in use, make sure the water tank is dry and inverted in the storage position.

WATER HOSE

- Wipe the water hose with a cloth.
- Do not bend, twist or pull the water hose. Otherwise the water hose may malfunction, e.g. stop absorbing water during operation.

WATER TANK PLUG

Wipe off dirt and dust from the plug with a cloth.

JET TIP

- Press the jet tip unlock button to remove the tip. Clean the tip with water and wipe with a cloth.
- It is recommended to replace the tip with a new one every 3 months, otherwise the effectiveness of the teeth cleaning will be reduced.

PROBLEM SOLVING

Problem	Cause	Solution
Jet tip cracking.	Serious damage to the jet tip.	Replace the tip with a new one matched to the model.
Water leaks from the water tank.	Seal is damaged.	Replace seal at an authorised service centre.
	Water remains in the water tank.	Empty the water tank after each use.
Leakage between the jet tip and the mounting hole.	The jet tip is not installed correctly.	Remove the jet tip and reinstall it.
No water flow or insufficient pressure.	The device is held in the wrong position.	Place the device in the correct position: the jet tip must be higher than the water tank.
	Insufficient water in the water tank.	Fill the water tank with water.
	Incorrect operating mode.	Select "Strong" mode.

INFORMATION FOR USERS ON HOW TO DEAL WITH WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT.

If the device, packaging, instruction manual and the like are marked with the symbol of a crossed-out wheeled waste container, it means that the device is subject to separate waste collection in accordance with Directive 2012/19/EU of the European Parliament and of the Council. The used device should not be disposed of with other household waste and should not be treated as municipal waste. By disposing of electro-waste in the container you pose a threat to the environment. Take the used device to a separate collection point for electrical and electronic equipment organized by the public administration. By segregating and submitting used electrical and electronic equipment for processing, recovery, recycling and disposal, you protect the environment from pollution and contamination, contribute to reducing the use of natural resources and lowering the cost of producing new equipment. Proper disposal and scrapping helps eliminate the adverse effects of scrapped products on the environment and human health. For details on the recyclability of this product, please contact your local municipality, city cleaning services or the store where you purchased the product.

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